

DEDICATED TO REVOLUTIONIZING
THE POWER CHAIR EXPERIENCE



BRYAN ANDERSON

Iraq War Veteran
and Quantum Rehab®
Spokesperson

QUANTUM® OUTBACK

FDA Class II Medical Device



quantumrehab.com



5	INTRO
6	iLEVEL®
8	OUTBACK
10	EDGE® 3
12	STRETTO
14	R-TRAK
16	4FRONT® 2
20	BARIATRIC LINE
22	TRU-BALANCE® 4
24	SEATING OPTIONS
26	STEALTH PRODUCTS® CUSHIONS
30	ELECTRONICS AND ACCESSORIES
34	COLOR OPTIONS
36	APPS
38	SERVICE AND SUPPORT
40	SPECIFICATIONS



Quantum® Power Chairs are
FDA Class II Medical Devices
designed to aid individuals with
mobility impairments.



FDA Class II
Medical Device



WELCOME TO QUANTUM®

Quantum Rehab® was born out of the desire to delight customers with the most advanced, consumer-inspired complex rehab power wheelchairs and related technologies possible.

Clinically-focused and consumer-inspired, Quantum® products are engineered and manufactured to meet the widest range of complex rehab needs with functional, durable, stylized, high-performance designs to increase daily functionality and quality of life. At Quantum, consumer needs and wishes are the driving force.

We're dedicated to not just meeting medical and clinical needs, but also quality-of-life needs. From the most advanced power seating for pressure management to USB ports, Bluetooth® and fender lights, no consumer need is overlooked.

Headquartered in Duryea, Pennsylvania, Quantum also has operations in Australia, Canada, China, Italy, France, Germany, the Netherlands, New Zealand, and the United Kingdom.

quantumrehab.com



Quantum® Power Chairs are
FDA Class II Medical Devices
designed to aid individuals with
mobility impairments.



FDA Class II
Medical Device



FDA Class II
Medical Device





"COOKING IS SO MUCH EASIER NOW THAT I CAN REACH THINGS. I CAN'T IMAGINE MY LIFE WITHOUT iLEVEL."

STEPHANIE WOODWARD
ROCHESTER, NY

FDA Class II
Medical Device

iLEVEL®

POWER SEAT ELEVATION

Standing Height + Speed + Stability = Authentic iLevel (like no other)

iLevel® power seat elevation technology features a revolutionary new power chair design for a better experience. iLevel improves daily functionality and quality of life by allowing operation of the power chair with the seat elevated while at walking speed.

FEATURES

- 12" of power adjustable seat height while at walking speed; maximum of 10" height on the Q6 Edge® HD
- iLevel's patented Extra Stability Technology® enhances safety for transfers, reach, and numerous other activities of daily living
- Available on the Stretto, Edge® 3 and Q6 Edge® 2.0
- Front and rear LED fender lights come standard with iLevel®



FDA Class II
Medical Device

Safer transfers



Improved reach



Social interaction at iLevel®



FDA Class II
Medical Device

FDA Class II
Medical Device

FDA Class II
Medical Device

MAKE YOUR OWN PATH.



"THE OUTBACK LETS ME LIVE LIFE THE WAY I WANT. WHETHER I'M WORKING ON THE FARM OR GOING FOR A HIKE WITH MY FRIENDS, THE OUTBACK ALLOWS ME TO GO WHEREVER I WANT TO GO."

BRYAN ANDERSON
NASHVILLE, TN

QUANTUM®
OUTBACK

4X4 ALL-TERRAIN

Quantum Rehab® welcomes adventurers and nature enthusiasts to unleash their wild side with the release of our cutting-edge 4x4 all-terrain power chair, the Quantum® Outback! Designed for those who dare to defy conventional limits, this rugged mobility solution promises to elevate any outdoor adventure, whether navigating wooded trails or exploring uncharted terrain.

SPECS


MAXIMUM SPEED: Up to 5 mph (6 mph with optional upgrade)

WEIGHT CAPACITY: 300 lbs. (350 lbs. with Q-Captain Seating)

TURNING RADIUS: 49.2"

POWER SEAT ELEVATION: Optional 10" safe seat elevation

COLORS

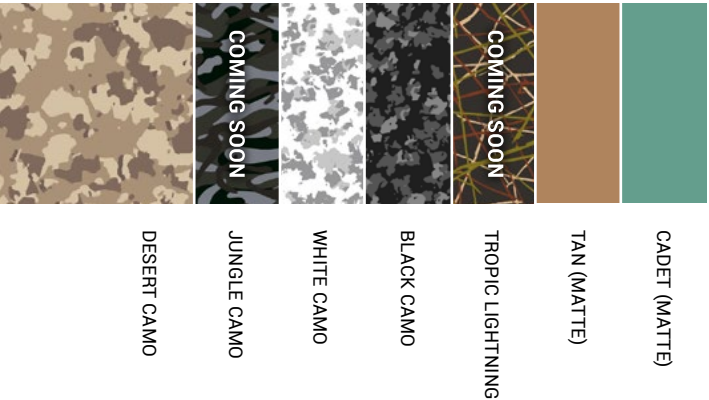

Available with TRU-Balance® 3 Power Positioning (shown) or Q-Captain Seating

Optional Quantum® backup camera

OPTIONAL
10" safe seat elevation at up to 3.2 mph

Rear LED lights

Q-Logic 3e Controller with Bluetooth®

Stealth Products® cushion featuring Coolcore® technology

Front LED lights

Independent SRS (Smooth Ride Suspension)

15" pneumatic tires for superb control over rough terrain




The Quantum® Outback is an FDA Class II Medical Device designed to aid individuals with mobility impairments.

FRONT LED LIGHTING

REAR LED LIGHTING

Q-LOGIC 3e CONTROLLER WITH BLUETOOTH®

INDEPENDENT SMOOTH RIDE SUSPENSION

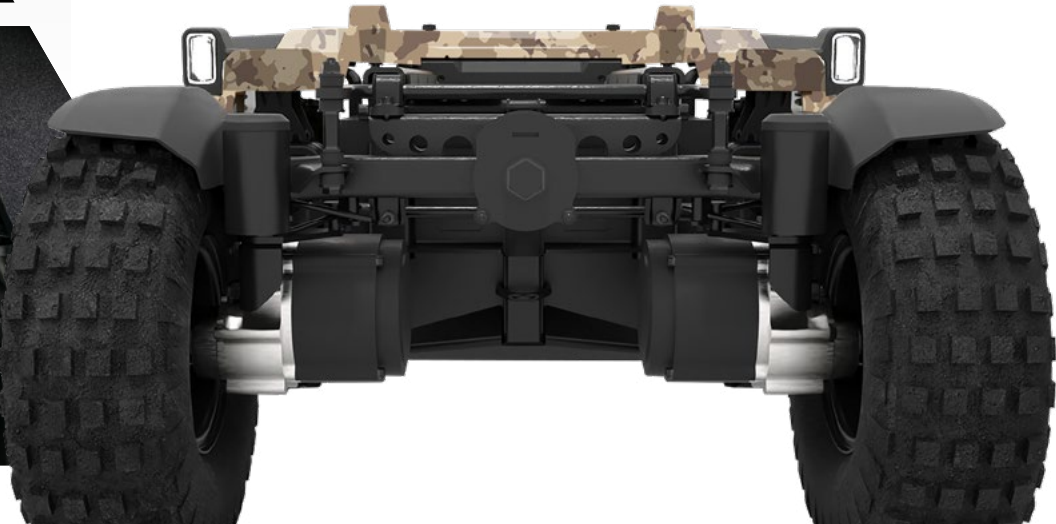

PROUDLY MADE IN
AMERICA
OF FOREIGN AND DOMESTIC PARTS



"THE EDGE 3 WITH iLEVEL HAS RESTORED MY INDEPENDENCE AT BOTH WORK AND HOME. iLEVEL HAS ALLOWED ME TO GET MORE INVOLVED AND START A DIALOGUE WITH PEOPLE, SOMETHING I WOULDN'T HAVE DONE BEFORE."

MERLISHA HENDERSON
PHOENIX, AZ

EDGE® 3

MID-WHEEL DRIVE

The Edge® 3 offers upgraded SRS (Smooth Ride Suspension), which provides a smooth, comfortable ride and advanced stability. Optional iLevel® is available at 4.5 mph with up to 12" of power adjustable seat height.

SPECS

MAXIMUM SPEED: 6.25 mph (7 mph with optional 4-pole motors and Q-Logic 3 / 4.5 mph when elevated with iLevel®)

WEIGHT CAPACITY: 300 lbs. (See Q6 Edge® HD for increased capacity)

TURNING RADIUS: 20.5"

POWER SEAT ELEVATION: Optional 12" iLevel

COLORS

* See website for a complete list of color options.



GREEN MACHINE
RED CORVETTE
DEEP PURPLE
LEMON CRUSH
RASPBERRY BERET
EVOLUTION ORANGE
BLUE JEAN

OPTIONAL
TRU-Balance® 4 Power Positioning Systems provide a choice of 10, 20, or 30 degrees of power anterior tilt.

OPTIONAL
12" iLevel® power seat elevation at 4.5 mph (standard on TRU-Balance® 4)

Featuring Extra Stability Technology® for safer transfers, reach and other activities of daily living

Q-Logic 3 Advanced Drive Controls with Bluetooth®

Stealth Products® cushion featuring exclusive Coolcore® technology

Upgraded SRS (Smooth Ride Suspension)



The Edge® 3 is an FDA Class II Medical Device designed to aid individuals with mobility impairments.

USB charger (standard)

Front and rear LED fender lights for seeing and being seen (standard)

Drive wheel color accents (standard)



PROUDLY MADE IN
AMERICA
OF FOREIGN AND DOMESTIC PARTS

NARROWEST POWER
BASE IN AMERICA AT
20.47"
WIDTH WITH 12.5"
DRIVE WHEELS



"TIGHT ROOMS, NARROW DOORWAYS, OUT-OF-REACH CABINETS AND VAN RAMPs... NOW I CAN GO WHERE I WANT TO GO, ALL BECAUSE OF THE STRETTO!"

BRYSON FOSTER
CONCORD, NC

OPTIONAL
TRU-Balance® 4 Power Positioning Systems provide a choice of 10, 20, or 30 degrees of power anterior tilt.

OPTIONAL
12" iLevel® power seat elevation at 3.5 mph (standard on TRU-Balance® 4)

Toolless removal of rear shroud with battery access door

Independent SRS (Smooth Ride Suspension)

Q-Logic 3 Advanced Drive Controls with Bluetooth® and USB charger

Stealth Products® cushion featuring exclusive Coolcore® technology

Front and rear LED fender lights



The Stretto is an FDA Class II Medical Device designed to aid individuals with mobility impairments.

STRETTO

MID-WHEEL DRIVE

The Edge 3 Stretto® is the narrowest, most maneuverable power base in America! In addition to its narrow width, it is equipped with independent SRS (Smooth Ride Suspension), making it ideal for children, teenagers and small adults. The Stretto is available with 12.5" or 14" drive wheels and optional 12" of iLevel® power adjustable seat height.

SPECS

MAXIMUM SPEED: 6.25 mph (3.5 mph when elevated with iLevel®)

WEIGHT CAPACITY: 300 lbs.

TURNING RADIUS: 21"

POWER SEAT ELEVATION: Optional 12" iLevel

COLORS

* See website for a complete list of color options.



BLUE JEAN

RED CORVETTE

DEEP PURPLE

LEMON CRUSH

RASPBERRY BERET

EVOLUTION ORANGE

GREEN MACHINE

Independent SRS (Smooth Ride Suspension)

The Stretto is the ideal pediatric power chair solution

Toolless access to electronics

THE EVOLUTION OF REAR WHEEL DRIVE



The R-TRAK is an FDA Class II Medical Device designed to aid individuals with mobility impairments.



R-TRAK

As a global innovator, we are creating a product that offers confidence and stability elevated at the speed of life. Our state-of-the-art facility has a robust team of experienced technicians who custom build each R-TRAK specifically to the individual consumer's needs.



Scan to learn about
the design and
development of
the R-TRAK.



R-TRAK

REAR-WHEEL DRIVE

The R-TRAK revolutionizes the rear-wheel drive power chair experience with independent SRS suspension on all four wheels, providing a smoother, more maneuverable ride. Standard features include fender lights, USB charger and drive wheel color accents. The R-TRAK is also available with our 10" safe seat elevation, allowing the client to drive up to 3.5 mph while elevated.

SPECS

MAXIMUM SPEED: Up to 7 mph (3.5 mph when elevated)


WEIGHT CAPACITY: 300 lbs.

TURNING RADIUS: 25.7"

POWER SEAT ELEVATION: Optional 10" safe seat elevation

COLORS

* See website for a complete list of color options.

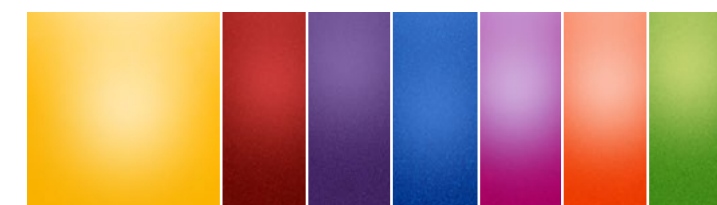

LEMON CRUSH

RED CORVETTE

DEEP PURPLE

BLUE JEAN

RASPBERRY BERET

EVOLUTION ORANGE

GREEN MACHINE

Independent SRS (Smooth Ride Suspension)

LED lighting in front caster
arms and rear fenders

Optional Quantum backup camera
can be mounted anywhere



PROUDLY MADE IN
AMERICA
OF FOREIGN AND DOMESTIC PARTS

450-LB.
WEIGHT CAPACITY

available on the
4FRONT® 2 HD

Optional backup camera with
3.5" display can be mounted
in multiple locations

OPTIONAL

TRU-Balance® 4 Power
Positioning Systems provide
a choice of 10°, 20°, or 30°
of power anterior tilt.

OPTIONAL

10" power adjustable safe
seat elevation at 3.5 mph
(standard on TB4)

Redesigned caster arms
for improved durability

STC (Smart
Traction Control)



Automotive-grade SRS
(Smooth Ride Suspension)

The 4Front® 2 is an FDA Class II Medical Device designed to aid individuals with mobility impairments.



"RIDING IN THE 4FRONT 2 IS UNLIKE ANY
CHAIR I HAVE EVER USED. THE 4FRONT 2
PUSHES THE LIMIT SET BY ITS PREDECESSOR
BEYOND MY EXPECTATIONS."

CHRISTIAN BUDNEY
SCRANTON, PA

Q-Logic 3e controller
featuring a magnetic
charging port (Q-Logic 3
also available)

Stealth Products®
cushion featuring
exclusive Coolcore®
technology

Front and rear LED
fender lights for
seeing and being
seen (standard)



4FRONT® 2

FRONT-WHEEL DRIVE

Featuring STC (Smart Traction Control)

The 4Front® 2 is the next generation of the 4Front® Power Chair! New Smart Traction Control (STC), along with Smooth Ride Suspension (SRS), provide stability for smooth obstacle transitions and unmatched driving performance. Additional new features include 4-pole motors and redesigned caster arms for improved durability and aesthetics and reduced weight.

S P E C S

MAXIMUM SPEED: 6.25 mph (3.5 mph when elevated)

WEIGHT CAPACITY: 300 lbs. (450 lbs. on 4Front® 2 HD)

TURNING RADIUS: 25.9"

POWER SEAT ELEVATION: Optional 10" safe seat elevation

COLORS

* See website for a complete list of color options.



EVOLUTION ORANGE

RED CORVETTE

DEEP PURPLE

LEMON CRUSH

RASPBERRY BERET

GREEN MACHINE

BLUE JEAN

450-lb. weight capacity is available
on the **4Front® 2 HD**

Automotive-grade SRS (Smooth Ride Suspension)

LED front and rear fender lights
(front: white / rear: red)



FDA Class II
Medical Device



Q6 EDGE® 2.0

MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 6.25 mph (3.5 mph when elevated)


WEIGHT CAPACITY: 300 lbs.

TURNING RADIUS: 20.5"

POWER SEAT ELEVATION: Optional 12" iLevel®

COLORS

* See website for a complete list.

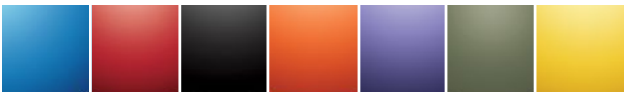

FDA Class II
Medical Device



AVAILABLE WITH
10" SAFE SEAT
ELEVATION AND
Q4 POWER TILT

J4

MID-WHEEL DRIVE

SPECS


MAXIMUM SPEED: 5 mph

WEIGHT CAPACITY: 300 lbs.

TURNING RADIUS: 21"

POWER SEAT ELEVATION: Optional 10" Safe Seat Elevation

COLORS

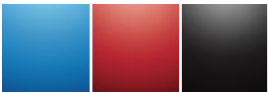

FDA Class II
Medical Device



Q6 EDGE® HD

BARIATRIC MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 4.7 mph (3.5 mph when elevated)

WEIGHT CAPACITY: 450 lbs. (400 lbs. with iLevel®)

TURNING RADIUS: 19.75"

POWER SEAT ELEVATION: Optional HD 10" iLevel®

COLORS

* See website for a complete list.



FDA Class II
Medical Device



Q1450

BARIATRIC FRONT-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 4.5 mph

WEIGHT CAPACITY: 600 lbs.

TURNING RADIUS: 26.5"

COLORS



J4E
shown with
optional 10"
Safe Seat
Elevation



FDA Class II
Medical Device

FDA Class II
Medical Device



AVAILABLE WITH
10" SAFE SEAT
ELEVATION AND
Q4 POWER TILT

J4E

MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 5 mph

WEIGHT CAPACITY: 300 lbs.

TURNING RADIUS: 20.5"

POWER SEAT ELEVATION: Optional 10" Safe Seat Elevation

COLORS



FDA Class II
Medical Device



FEATURING
Q4 HEAVY DUTY
POWER TILT

J4 HD

MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 5 mph

WEIGHT CAPACITY: 400 lbs.

TURNING RADIUS: 21"

COLORS



FDA Class II
Medical Device



J6

MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 4 mph

WEIGHT CAPACITY: 300 lbs.

TURNING RADIUS: 22"

COLORS



#1 FOR REHAB POWER

BUILT TOUGH

Quantum® bariatric power chairs all offer an innovative design with superior outdoor performance and superb indoor maneuverability. Our baratric line of power chairs offers speeds as high as 5 mph and weight capacities up to 600 lbs.

All of our bariatric power chairs can be equipped with power positioning for greater comfort, such as TRU-Balance® HD Power Tilt. This seating system can assist with pressure relief. To provide additional support, Stealth® cushions and backs are available to provide lumbar support or coccyx and sacral relief. All Quantum bariatric power chairs are available with Q-Logic environmental controls featuring Bluetooth® technology.



Quantum® Power Chairs are FDA Class II Medical Devices designed to aid individuals with mobility impairments.

	Q6 Edge® HD	Q1450	J4 HD	4Front® 2 HD
Capacity	450 lbs.	600 lbs.	400 lbs.	450 lbs.
Top Speed	4.7 mph	4.5 mph	5 mph	5 mph
Turning Radius	19.75"	26.5"	21"	24.5"
Power Elevation	Yes (10")	N/A	N/A	Yes (10")

** Please consult the back of the catalog for detailed specifications.



TRU-Balance® 4 Power Positioning Systems

WITH ANTERIOR TILT AND MEMORY SEATING

Maximum weight capacity of up to 300 lbs.

Up to 8 customizable memory seating positions recall the saved positions of all actuators, offering the choice between sequential or synced movement for each seat function.

Memory seating profiles can be labeled for easy recall of most used positions.

Seat elevation combined with anterior tilt provides optimal access.

Consumers may find it easier to perform ADLs such as: accessing sinks and cooktop surfaces, meal preparation and reaching items in cabinets.

TRU-Balance® 4 with Dual Actuator AFP shown on the Edge 3 Stretto®



DUAL ACTUATOR POWER ARTICULATING FOOT PLATFORM

Lowers to assist with transfers

Articulation range of 0-7.25"

Elevation range of 85-165 degrees

Only available on TRU-Balance® 4



FDA Class II Medical Device

Chest Bar



Knee Blocks



TRU-BALANCE® 4

POWER POSITIONING SYSTEMS

FEATURES

- Anterior tilt options: 10, 20 and 30 degrees
- Up to eight customizable memory seating positions
- Chest bar and knee blocks available (knee blocks and chest support required with 20 or 30 degrees of anterior tilt)
- Available on the Edge 3 Stretto®, 4Front® 2, and Edge® 3
- Optional Accu-Trac Advanced Tracking Technology
- Optional WC19 crash-test approved occupied transit kit option
- Available seat sizes: 16-22" width and 14-22" depth

POWER RECLINE

- 85° - 165° of recline
- 300 lbs. weight capacity
- Available with power tilt and power adjustable seat lift
- No increase in seat-to-floor height
- Standard shear reduction and cantilever armrests

POWER POSTERIOR TILT

- 0° - 50° of tilt
- Lower seat-to-floor height
- Unmatched adjustment and durability
- Maintains center of gravity shift

DRIVE SPEEDS IN ANTERIOR TILT

- Up to 10° of anterior tilt, up to full lift: **2 MPH**
- 10-20° of anterior tilt, up to 6" of lift: **0.5 MPH**
- 10-20° of anterior tilt, over 6" of lift: **DRIVE LOCKOUT**
- Over 20° of anterior tilt: **DRIVE LOCKOUT**

COLORS

* See website for a complete list of color options.



LEMON CRUSH

RED CORVETTE

DEEP PURPLE

EVOLUTION ORANGE

RASPBERRY BERET

GREEN MACHINE

BLUE JEAN

TRU-Balance® 4 Power Recline



TRU-Balance® 4 Power Posterior Tilt



TRU-Balance® 4 with 12" iLevel® Power Elevation





"I'M ABLE TO ELEVATE TO EVERYONE'S HEIGHT WITH iLEVEL. I FEEL MORE INCLUDED."

MORGAN STEWARD
COVINGTON, GA

Morgan Steward shown on his Stretto with TRU-Balance® 3 power positioning.

FDA Class II
Medical Device

SEATING AND POSITIONING



TRU-BALANCE® FLEX SEATING

No-charge static seat for the Edge 3 Stretto®, Edge® 3, Q6 Edge® 2.0, J4 and Rival®

Easily adjustable, versatile design that can grow with the user

Same ease of width and depth adjustment as TRU-Balance® 3

Seat sizes 15"-22" wide and 15"-20" deep



CONTOURED SEATING

Black vinyl captain's seat with optional armrests

Compatible with TRU-Balance 3 power articulating foot platform

Available in high back, medium back, solid seat pan, limited recline, and depth-adjustable



PEDIATRIC SEATING

Easily adjustable, this versatile design can grow with the user

Several width/depth adjustment packages start at 10"

Accepts a wide variety of after-market seats and backs

Available in numerous shroud colors that kids will love

TRU-Balance®
Flex Seating



Contoured Seating



Pediatric Seating





FEATURING EXCLUSIVE
COOLCORE
CHEMICAL-FREE COOLING TECHNOLOGY



30+ UPF



Dry



Chemical
Free



Breathable
Comfort

"While the other makers of cooling fabrics use a number of chemicals to enhance wicking, transportation and evaporation, Coolcore fabrics are chemical free. To make that possible, we changed the essence of cooling fabric technology right down to its very core."

COOLCORE®


STEALTH PRODUCTS®

CUSHIONS

Stealth Products® Cushions feature patented Coolcore® technology, a chemical free, four-way stretch material that offers moisture wicking, moisture transportation, and moisture evaporation properties. Our cushions also use a color guide that identifies the suggested HCPC code and function family by a unique patch.

FEATURES


- Chemical free
- Regulates temperature
- Manages and distributes your body heat
- Ships out within 2 days on or off chair

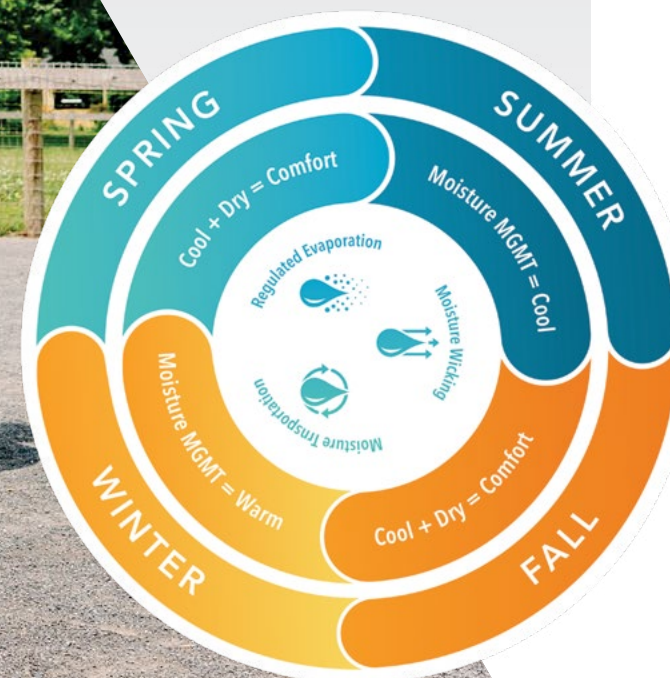

FDA Class II
Medical Device



Convenient built-in side/front pouches

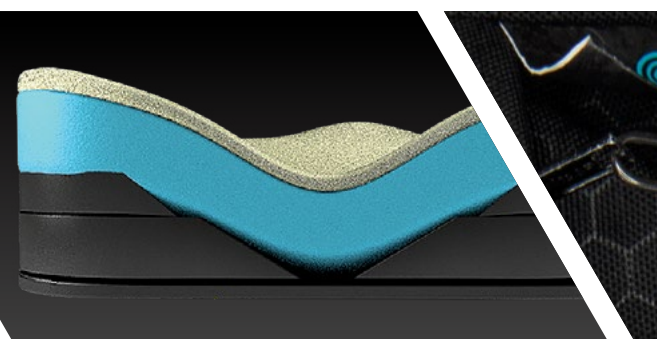









Constructed with FRG Sunmate
flame-retardant material



Easy-pull rubber loop
zipper handles





	TRU-Comfort 2 SPP	Essence SPP	Solution® SPP	Spectrum® Foam/Gel SPP	Premiere P	Glacial SP	Zen SP	Simplicity G
DESCRIPTION ¹	 <p>Multi-layered foam</p> <p>Deep contour, CoolCore® cushion cover has “pockets” on sides</p>	 <p>Low profile Tru-Comfort 2, lateral and pre-ischial shelf</p> <p>Contoured, HD molded foam base, coccyx and sacral relief</p>	 <p>Standard viscoelastic foam top with twin cell aqueous viscous gel (heavier gel) insert in ischial well, mid-level contour</p>	 <p>SunMate FRG foam top (SunMate FRG XXSoft (blue); Medsoft (green)) with multi-layered foam insert or twin cell ultra-lightweight viscous gel insert in ischial well, sacral relief</p>	 <p>Gray X-Soft SunMate FRG foam topper, high density positioning wedges built in for offloading.</p>	 <p>New SIM G base, fluid resistant overlay, top and sides</p> <p>TPU heat-sealed, low profile, lightweight cooling gel topper</p>	 <p>New SIM G base, SunMate FRG XXSoft (blue topper), low profile, lightweight</p>	 <p>High-density foam</p> <p>Stealth-Tek cover standard</p> <p>More contour to sacral/coccyx relief well</p>
WEIGHT CAPACITY	14-15” W: 175 lbs. 16-21” W: 350 lbs. 22+” W: 450 lbs.	10-15” W: 175 lbs. 16-21” W: 350 lbs. 22” W: 450 lbs.	10-15” W: 175 lbs. 16-20” W: 350 lbs. 22+” W: 500 lbs.	14” W: 175 lbs. Foam & Gel 16-20” W: 350 lbs. Foam & Gel 22+” W: 650 lbs. Foam; 550 lbs. Gel	14-15” W: 175 lbs. 16-21” W: 350 lbs. 22” W: 500 lbs.	10-15” W: 175 lbs. 16-21” W: 350 lbs. 22+” W: 400 lbs.	10-15” W: 175 lbs. 16-21” W: 350 lbs. 22+”: 400 lbs.	10-15” W: 175 lbs. 16-21” W: 350 lbs. 22+” W: 400 lbs.
OVERALL WEIGHT	4 lbs.	3.5 - 5 lbs.	5.4 lbs.	Foam: 3 lbs. Gel: 3.5 lbs.	3.5 lbs.	3.5 lbs.	2.5 lbs.	2 lbs.
HEIGHT	4.5”	3.5”	3”	3”	3.5”	2.5”	3”	2”
DIMENSIONS ²	14x14” up to 24x24” Odd sizes in 15-21” W SOF available	10x10” up to 22x24” Odd sizes in 15-21” W SOF available	10x10” up to 24x24” Odd sizes in 15-19” W SOF available	Foam: 14x18” up to 32x28” Gel: 14x18” up to 24x26” Odd sizes in 17” and 19” W SOF available	14x13” up to 22x22” Odd sizes up to 21” W SOF available	10x10” up to 24x24” Odd sizes in 15-21” W No SOF	10x10” up to 24x24” Odd sizes in 15-21” W SOF available	10x10” up to 24x24” Odd sizes in 15-21” W No SOF
MSRP	\$420 - \$525	\$420 - \$625	\$394 - \$452	\$483 - \$557 (\$620 for Spectrum Foam SPP larger than 24”)	\$341 - \$499	\$347 - \$462	\$331 - \$446	\$116 - \$173
HCPCS	E2607/E2608	E2607/E2608	E2607/E2608	E2607/E2608	E2605/E2606	E2603/E2604	E2603/E2604	E2601/E2602
ALLOWANCES ³	\$332.23 / \$398.98	\$332.23 / \$398.98	\$332.23 / \$398.98	\$332.23 / \$398.98	\$308.53 / \$481.33	\$173.74 / \$215.97	\$173.74 / \$215.97	\$70.10 / \$136.85

1) CoolCore® is standard on all cushions except the Simplicity G.
 2) Special order sizes available for all cushions except the Glacial SP and Simplicity G. Contact sales for additional assistance.
 3) Allowances listed are estimations for cushions billed in conjunction with Group 3 power chair bases using the KU modifier and may vary by state.

DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement. HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. Payer coding, coverage, and bundling guidelines may apply.



Stealth® ADI Q-Back 2.0

FEATURES

- 450-lb. weight capacity
- CoolCore® Cover with 1.25” of Medium Soft Sunmate® FRG foam
- Two styles: AL Style (Aluminum) and AFT Style (Aluminum Flat Top)
- Available configurations: Active Contour (2.5” depth), Deep Contour (4” depth), and Extra Deep Contour (6” depth)
- Widths: 14”-22”
- Shell heights: Standard (16”), Tall (18”), Extra Tall (20”), Extra Extra Tall (22”)





"MY EDGE 3 GIVES ME THE OPPORTUNITY TO USE ADVANCED TECHNOLOGY THAT GIVES ME SO MUCH INDEPENDENCE."

ZOE HERNANDEZ



Q-LOGIC 3

ADVANCED DRIVE CONTROLS

The Q-Logic 3 Advanced Drive Controls System from Quantum Rehab® is like no other power wheelchair electronics system. Designed from clinician and consumer feedback, Q-Logic 3 is user-friendly and highly customizable.

FEATURES

- Wireless Bluetooth® programming through Windows devices
- Bluetooth access to computer, tablet, and phone functions
- iAccess programmable push button and toggle options for customized seating access
- Clinic Mode:
 - Quickly set up multiple input devices in client evaluations
 - Automatically assigned drive modes to input devices

The Q-Logic 3e features a magnetic charging port for clients with limited dexterity



Q-LOGIC 3 EX
ENHANCED DISPLAY




Q-LOGIC 3 EX
JOYSTICK CONTROLLER



Q-LOGIC 3e
JOYSTICK
CONTROLLER



Q-Logic 3
home screen

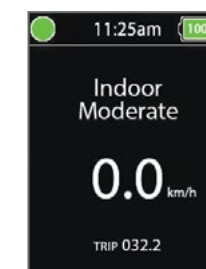


Q-Logic 3
drive screen

ENHANCED DISPLAY

The optional Q-Logic 3 EX Enhanced Display serves as a separate screen used with specialty drive controls to display system data. With the Enhanced Display's built-in environmental controls, Q-Logic 3 can operate a wide range of appliances that use infrared (IR) remote controls, such as televisions, and serves as a Bluetooth® computer mouse emulator.



Q-Logic 3 iAccess Programmable Switchbox
available with push buttons or toggles

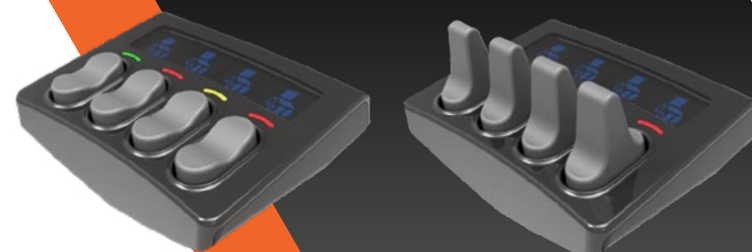



































Q-Logic 3 Stand
Alone Joystick



Q-Logic 3 Specialty Controls
Interface Module



NON-EXPANDABLE		JOYSTICK OPTIONS		SIP-N-PUFF		MULTI-SWITCH				
										
Q-Logic NE and NE+	Q-Logic NES			Stealth® i-Drive® Sip-N-Puff		Stealth® i-Drive® 4 Stealth Pro, 3-Switch Head Array	Stealth® i-Drive® Tri-Array, 3-Switch Head Array	Stealth® i-Drive® Curved Array		
EXPANDABLE (Q-LOGIC 3)				SINGLE SWITCH				TRAY ARRAYS		
										
Q-Logic 3 EX Joystick	Q-Logic 3e Joystick	Q-Logic 3 EX Enhanced Display	Q-Logic 3 Stand Alone Joystick	Q-Logic 3 iAccess Programmable Switchbox available with push buttons or toggles	Q-Logic 3 Specialty Interface Module	Stealth® Egg Switch	Mo-Vis™ Twister Basic Button	Mo-Vis™ Twister Satellite Button	Stealth® i-Drive® Eclipse Tray with Four Capacitive Proximity Sensors	Stealth® i-Drive® Eclipse Tray with Four Fiber Optic Proximity Sensors
PROPORTIONAL DRIVE CONTROLS						MISC. ACCESSORIES				
										
Mo-Vis™ All-Round Joystick	Mo-Vis™ All-Round Light Joystick	Mo-Vis™ Multi Joystick	Mo-Vis™ Foot Control		Stealth® i-Drive® Control Harness	Front and Rear LED Fender Lights	Hydration System	Transfer Bars (Pair)	Backup Camera	
										
Stealth® Precision Mini Proportional Joystick	Stealth® Mushroom Joystick	Mo-Vis™ Micro Joystick				Personal Item Hook	Phone Holder	Personal Items Pouch	Glove Box	

COLOR OPTIONS



BLUE JEAN (MATTE)
PPW-BLU-11
Available on Stretto, Edge® 3, 4Front® 2, 4Front® 2 HD and R-TRAK



RED CORVETTE
PPW-RED-01
Available on all Quantum® bases



MIDNIGHT BLUE
PPW-BLU-10
Available on Edge® 3/2.0/2.0X/Z/HD, Stretto, 4Front® 2, 4Front® 2 HD and R-TRAK



EVOLUTION ORANGE
PPW-ORN-07
Available on all Quantum® bases except for J4 and J4 HD



LEMON CRUSH
PPW-YEL-04
Available on Stretto, Edge® 3/2.0/2.0X, 4Front® 2, 4Front® 2 HD and R-TRAK



GREEN MACHINE
PPW-GRN-11
Available on Edge® 3/2.0/2.0X/Z/HD, Stretto, 4Front® 2, 4Front® 2 HD and R-TRAK



RASPBERRY BERET
PPW-PNK-04
Available on Stretto, Edge® 3, 4Front® 2, 4Front® 2 HD and R-TRAK



PRETTY IN PINK
PPW-PNK-03
Available on all Quantum® bases except for J4 and J4 HD



DEEP PURPLE
PPW-PRP-05
Available on, Edge® 3/2.0/2.0X/Z/HD, Stretto, 4Front® 2, 4Front® 2 HD and R-TRAK



ELECTRIC BLUE
PPW-BLU-09
Available on all Quantum® bases



NIGHT TRAIN (MATTE)
PPW-BLK-03
Available on Edge® 3/2.0/2.0X/Z/HD, Stretto, 4Front® 2, 4Front® 2 HD and R-TRAK



MILITARY GREEN (MATTE)
PPW-GRN-08
Available on Stretto, Edge® 3/2.0/2.0X, 4Front® 2, 4Front® 2 HD and R-TRAK



GREY STREET
PPW-SIL-03
Available on Edge® 3/2.0/2.0X/Z/HD, Stretto, 4Front® 2, 4Front® 2 HD and R-TRAK



BACK IN BLACK
PPW-BLK-01
Available on all Quantum® bases




MILITARY DIGITAL CAMO
HYDROGRAPHICS
Available on Stretto, Edge® 3, 4Front® 2 and 4Front® 2 HD



FADED GLORY
HYDROGRAPHICS
Available on Stretto, Edge® 3, 4Front® 2 and 4Front® 2 HD

QUANTUM® OUTBACK COLOR OPTIONS

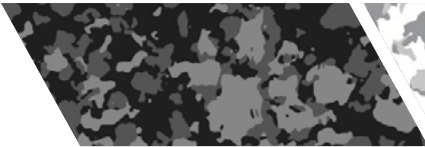

BLACK CAMO
DGN2301319-WCM9-PCL2



WHITE CAMO
DGN2301319-WCM10-PCL2



DESERT CAMO
DGN2301319-WCM11-PCL2



JUNGLE CAMO
DGN2301319-WCM13-PCL2
COMING SOON!




TROPIC LIGHTNING CAMO
DGN2301319-WCM12-PCL2
COMING SOON!



CADET
DGN2301319-TGY3-PCL2



TAN
DGN2301319-TTN1-PCL2

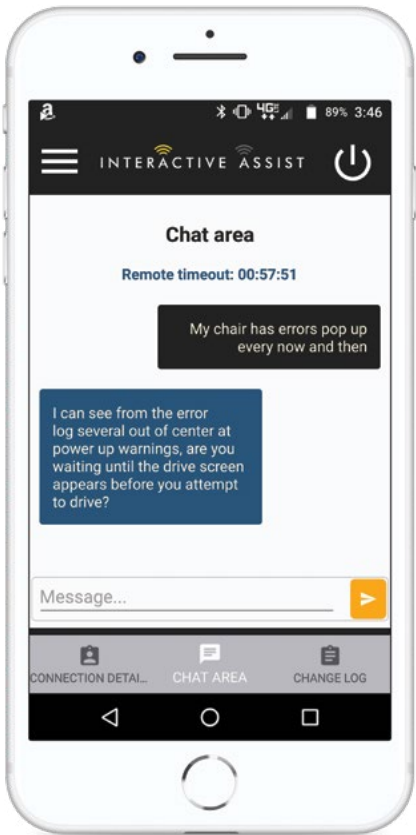

INTERACTIVE ASSIST

Reduce service visits with this innovative app.

- SCREEN MIRRORING
- REAL-TIME MONITORING
- EFFICIENCY IMPROVEMENTS
- TROUBLESHOOTING
- BATTERY MANAGEMENT
- DIAGNOSTIC HISTORY
- BASIC PROGRAMMING
- REMOTE ACCESS
- CHANGE LOG HISTORY
- SAVING XPA FILES
- MOTOR USAGE HISTORY
- TEXT CHATTING
- PROGRAMMING IR MENU
- ACCESS TO BLUETOOTH FEATURES
- DEALER-TO-DEALER PROGRAMMING ACCESS

Q-Logic 3 now features Interactive Assist. Interactive Assist from Quantum Rehab® is an app which provides a direct real-time connection from a power chair's electronics to a provider's programming station. The connection provides the technician with full system and diagnostics information along with tools, such as a real-time mirror image of the power chair's electronics display, to eliminate wasteful trips.



Download Interactive Assist today!



ECON-I

Convenience at your fingertips!

The Econ-I app is now available for download on Apple® iPhone® devices, allowing you to program Quantum® Power Chairs conveniently from your iPhone or iPad. Compatible devices include the iPhone 6S/6s Plus/SE or higher, and the iPad Air/iPad Mini 2/iPad Pro/iPad 5th generation or higher. All devices must support iOS version 12 or higher.



Download ECON-I today!

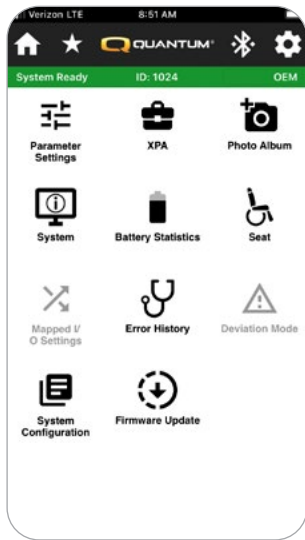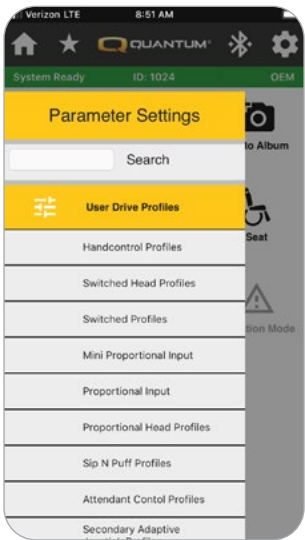


PQ-365

All Quantum® and Pride® product information combined into one app!




LIVE CHAT

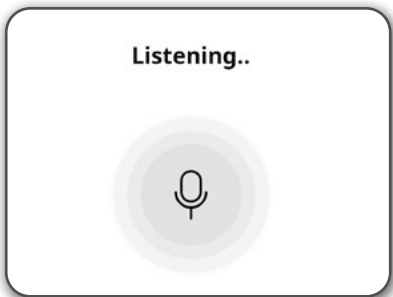


VOICE SEARCH



ZOOM FEATURE FOR PRODUCT IMAGES

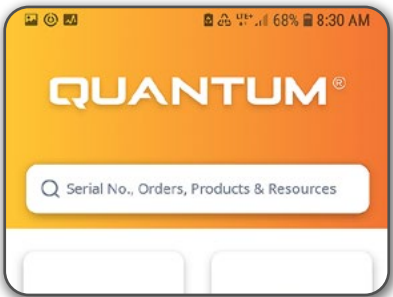

SERIAL NUMBER AND ORDER SEARCH

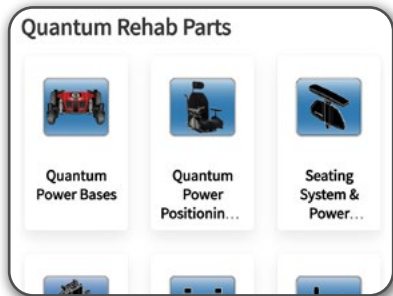

ILLUSTRATED PARTS BREAKDOWNS

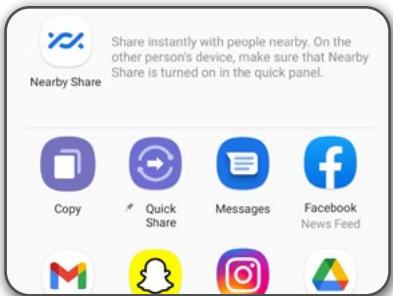


IMAGE SEND OPTION



Download PQ-365 today!

SERVICE AND SUPPORT

SALES SUPPORT 866.800.2002

Quantum Sales Support assists with design, ordering, pricing and quality assurance.

CUSTOM GROUP 877.465.2999

The Quantum Custom Group provides innovative custom solutions and specialized product modifications.

TECHNICAL SUPPORT EXT. 5200

Quantum Technical Support offers assistance with troubleshooting and parts order placement.

PRODUCT PLANNING & REIMBURSEMENT

Assistance with understanding Medicare codes and coverage criteria, product planning, provision guidelines, claims coding, pre-screens and assistance in regard to high-end rehab clients.


SET TO SPEC SERVICE

Set to Spec Service from Quantum Rehab® is designed to assist you in delivering Q6 Series Power Chairs to your clients. We will pre-set, to your specifications, seat-to-floor heights, seating center of gravity, armrest heights, front riggings, thoracic and/or thigh laterals.







LIVE LENS

Connect with technical service for remote, hands-on troubleshooting from anywhere in the field. Contact our Tech Service Dept. for details.

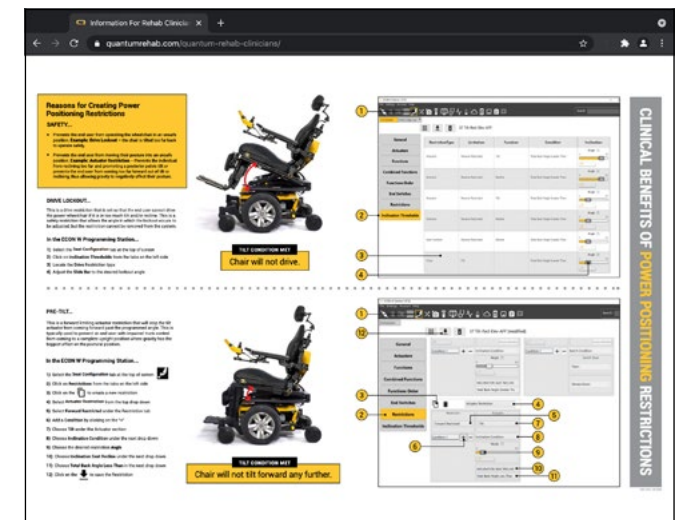


CLINICAL RESOURCES

The Clinical Resources section of the Quantum website is specifically tailored to clinicians' needs, featuring order forms, clinical references, and reimbursement information.

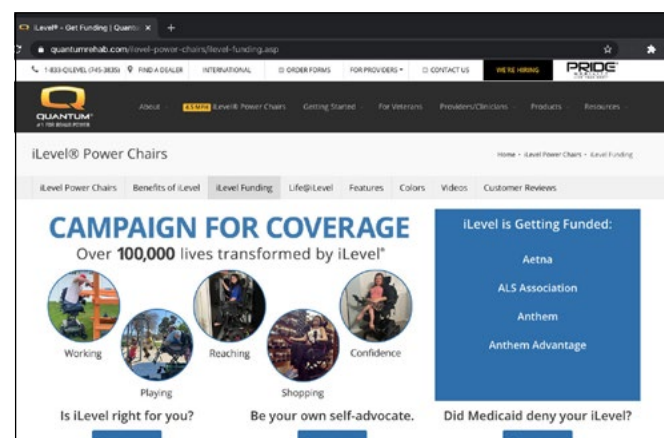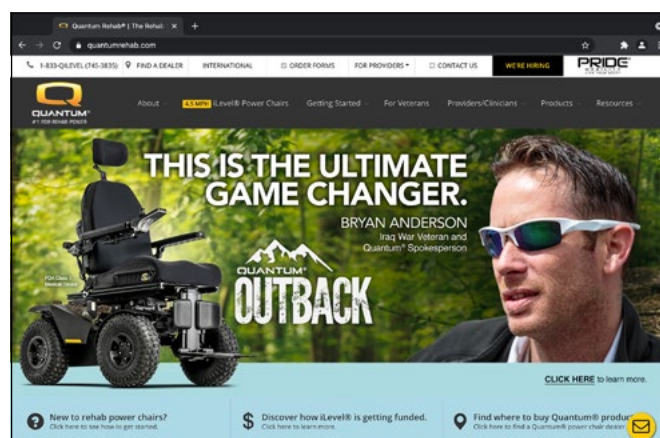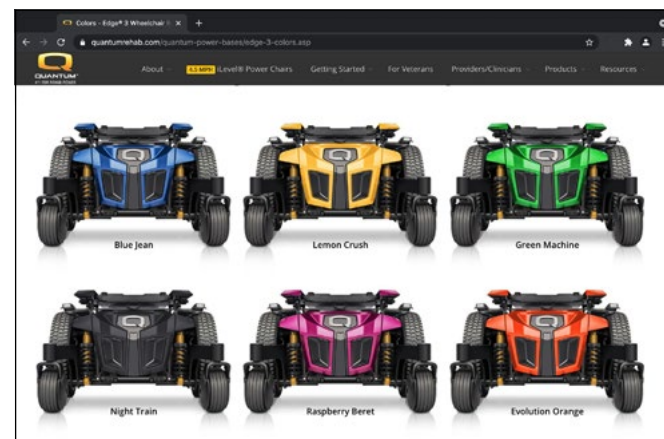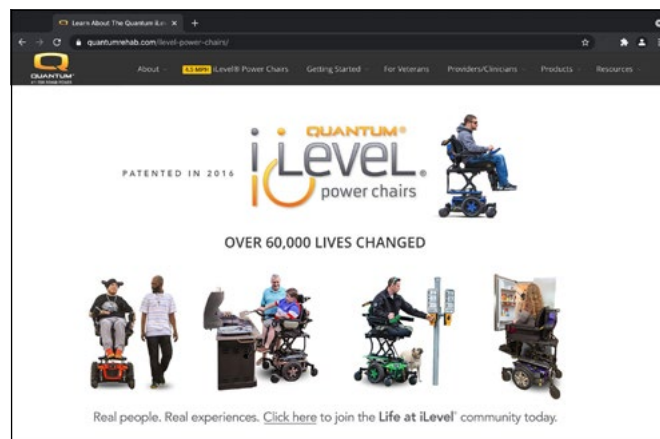

QUANTUMREHAB.COM

Quantum's official website is the ultimate source for everything you need to know about Quantum® products. Highlights include detailed product specifications, order forms, videos, user testimonials and much more.

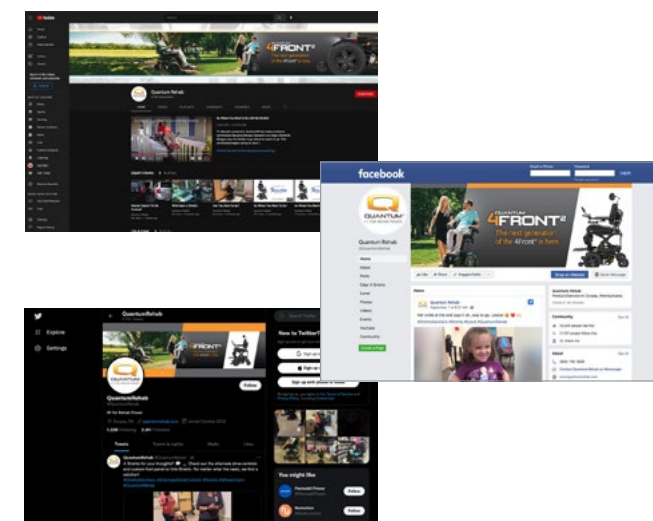

SOCIAL MEDIA

Stay up to date on the latest Quantum® events and product releases through Twitter, Facebook, YouTube, Instagram and TikTok.



STEALTHPRODUCTS.COM

Stealth Products® launched in 1999, in response to a need for versatile, high-quality head and neck positioning systems, as well as switch site location hardware. Their products include: head and body support/positioning, joystick mounts, trays, seating, cushions, and pediatric mobility products.

- 1) Standard

2) Optional

3) Requires HD Suspension

4) Speed varies with user weight, terrain type
battery charge, battery condition and tire
pressure.

5) Without front riggings

6) Battery capacity may change when adding
power positioning systems.

7) Weight capacity is dictated by the selection
of seating.

8) All measurements +/- 1/8" with 18" x 18" seat,
fixed towers, 300 lbs. weight capacity and solid
tires (add 3/8" with pneumatic tires). Depths
may increase seat-to-floor height. Contact
Customer Service for details. Listed
measurements are in the lowest position;
does not include elevating system unless
specified.

9) 250 lbs. weight capacity with TRU-Balance®

10) A minimum seating size of 20" x 20" is required
with iLevel® on the Q6 Edge® HD; maximum
seating size is 24" x 24"

11) Without foot riggings

12) All Quantum units can be ordered with
specialty controls

13) 6.25 mph speed available only with Q-Logic 3.

14) Battery weight may vary +/- 3 lbs. based on
manufacturer.

15) Reduced to 400 lbs. with iLevel® power
adjustable seat height

16) Seat-to-floor heights are measured with the
unit unloaded. Due to the 4Front's suspension
system, seat-to-floor heights are approximately
0.25" lower when the unit is loaded.

17) 7 mph can only be achieved with optional
4-pole motors and Q-Logic 3. Maximum
speed without 4-pole motors is 6.25 mph
with Q-Logic 3.

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement. The information contained herein is correct at the time of publication. We reserve the right to alter specifications without prior notice.

Quantum® Power Chairs are FDA Class II Medical Devices designed to aid individuals with mobility impairments.



	OUTBACK	STRETTO	EDGE® 3	EDGE® 2.0	EDGE® HD	J4
HCPCS CODES For more information on HCPCS code coverage, please consult our Product Planning & Reimbursement Center.	K0899 - Power Mobility Device, Not coded by DME PDAC	Group 3: Standard (K0848/K0849), Single Power (K0856/K0857), Multiple Power (K0861)	Group 3: Standard (K0848/K0849), Single Power (K0856/K0857), Multiple Power (K0861)	Group 3: Standard (K0848/K0849), Single Power (K0856/K0857), Multiple Power (K0861)	Group 3: Heavy Duty (K0850/K0851), Single Power HD (K0858), Multiple Power HD (K0862)	Group 2: Standard (K0822/K0823), Single Power (K0835)
AVAILABLE ON VA CONTRACT	Coming Soon to FSS Contract	FSS Contract V797D-70001	FSS Contract V797D-70001	FSS Contract V797D-70001	FSS Contract V797D-70001	FSS Contract V797D-70001
WEIGHT CAPACITY ⁷	300 lbs. (350 lbs. w/ Q-Captain)	300 lbs.	300 lbs.	300 lbs.	450 lbs. ¹⁵	300 lbs.
MAXIMUM SPEED ⁴	Up to 5 mph (6 mph with high speed package)	Up to 6.25 mph (3.5 mph with iLevel®) ¹³	Up to 6.25 mph ¹³ Optional: Up to 7 mph ¹⁷ (4.5 mph with iLevel®)	Up to 6.25 mph ¹³ (3.5 mph with iLevel®)	Up to 4.7 mph (3.5 mph with iLevel®)	Up to 5 mph
GROUND CLEARANCE	4.92" (at motors)	2.9"	2.7" (battery tray)	2.7" (battery tray)	3" (base frame)	2.9" (battery tray)
TURNING RADIUS ⁵	49.2"	21"	20.5"	20.5"	19.75"	21"
LENGTH ¹¹	42"	36"	35.5"	35.5"	35.6"	35.5"
WIDTH	28.82"	21.5" (14" drive wheels) 20.47" (12.5" drive wheels)	24.4"	24.4"	26.25"	24.5"
DRIVE WHEELS	15" pneumatic	14" or 12.5" pneumatic or flat free tires; Optional: 14" knobby	14" Solid Optional: 14" Pneumatic; 14" Knobby	14" Solid Optional: 14" Pneumatic; 14" Knobby	14" Run Flat Optional: 14" Pneumatic; 14" Knobby	14" Flat Free
CASTER WHEELS	N/A	Front: 5" flat free Rear: 6" flat free	Front: 6" Solid Rear: 6" Solid	Front: 6" Solid Rear: 6" Solid	Front: 6" Solid Rear: 6" Solid	Front: 6" Solid Rear: 6" Solid
ANTI-TIP WHEELS	N/A	N/A	N/A	N/A	N/A	N/A
BASE WEIGHT	168.6 lbs.	126.4 lbs.	138.6 lbs.	138.6 lbs.	143.5 lbs.	131.6 lbs.
BATTERY WEIGHT ¹⁴	58.6 lbs.	40 AH: 31 lbs. each NF-22: 40 lbs. each	40 lbs. each	40 lbs. each	52 lbs. each	40 lbs. each
DRIVE WHEEL CONFIGURATION	4 Motor, Independent Drive	Mid-wheel Drive	Mid-wheel Drive	Mid-wheel Drive	Mid-wheel Drive	Mid-wheel Drive
MOTOR PACKAGES	4-pole	2-pole	2-pole (Optional: 4-pole)	2-pole	4-pole	2-pole
BATTERIES ⁶	Group 27	40AH (2) or NF-22 (2)	NF-22 (2)	NF-22 (2)	Group 24 (2)	NF-22 (2)
BATTERY CHARGER	8A Off-board	8A Off-board	8A Off-board	8A Off-board	8A Off-board	8A Off-board
AVAILABLE ELECTRONICS ¹²	120A Q-Logic 3	120A Q-Logic 3; 75A Q-Logic NE; 75A Q-Logic NE+	120A Q-Logic 3; 75A Q-Logic NE; 75A Q-Logic NE+	120A Q-Logic 3; 75A Q-Logic NE; 75A Q-Logic NE+	120A Q-Logic 3; 75A Q-Logic NE; 75A Q-Logic NE+	75A Q-Logic NE; 75A Q-Logic NE+
SEATING SIZES						
TRU-BALANCE® 4	N/A	W: 16-22" D: 14-22"	W: 16-22" D: 14-22"	N/A	N/A	N/A
SYNERGY®	N/A	W: 10-22" D: 10-22"	W: 10-22" D: 10-22"	W: 10-22" D: 10-22"	W: 20-28" D: 15-28"	N/A
TRU-BALANCE® 3 STATIC	N/A	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	N/A	N/A
TRU-BALANCE® 3 POWER TILT	W: 16-22" D: 16-20"	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	N/A	N/A
TRU-BALANCE® 3 LIFT & TILT	W: 16-22" D: 16-20"	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	N/A	N/A
TRU-BALANCE® 3 HD STATIC	N/A	N/A	N/A	N/A	W: 20-32" D: 20-25"	N/A
TRU-BALANCE® 3 HD POWER TILT	N/A	N/A	N/A	N/A	W: 20-32" D: 20-25"	N/A
TRU-BALANCE® 3 HD LIFT & TILT	N/A	N/A	N/A	N/A	W: 20-24" D: 20-25" ¹⁰	N/A
TRU-BALANCE® 2 HD STATIC	N/A	N/A	N/A	N/A	N/A	N/A
TRU-BALANCE® 2 HD POWER TILT	N/A	N/A	N/A	N/A	N/A	N/A
TRU-BALANCE® FLEX SEAT	N/A	W: 15-22" D: 15-20"	W: 15-22" D: 15-20"	W: 15-22" D: 15-20"	N/A	W: 15-22" D: 15-20"
Q4 POWER TILT	W: 16-22" D: 16-22"	N/A	N/A	N/A	N/A	W: 16-20" D: 16-22"
SEAT-TO-FLOOR HEIGHTS ^{13,16}						
TRU-BALANCE® 4 ⁸	N/A	18-19"	17.5-19.25"	N/A	N/A	N/A
SYNERGY® / STATIC ⁸	N/A	16"-18" (Static only)	Synergy®: 15.875 - 20.875" Static: 16.5 - 18.5"	Synergy®: 15.875 - 20.875" Static: 16.5 - 18.5"	16.5-18.5"	N/A
POWER TILT ⁸	N/A	17.5-19.5"	17.5-19.5"	17.5-19.5"	18.5-20.5"	N/A
LIFT & TILT ⁸	N/A	18-20"	17.5-19.5"	17.5-19.5"	18.5-20.5"	N/A
TRU-BALANCE® FLEX SEAT ⁸	N/A	16.5-18.5"	16.5-18.25"	16.5-18.25"	N/A	16.5-18.25"
Q4 POWER TILT	N/A	N/A	N/A	N/A	N/A	17.25-19.25"
Q4 LIFT & TILT	N/A	N/A	N/A	N/A	N/A	N/A
Q CAPTAIN SEAT	Captain: 19.5-21.5"	SSP: 16.5-18.5" (Flex); 18-20" (Lift) Captain: 20.5-22.5" (Flex); 22-24" (Lift)	SSP: 16.5-18.25" (Flex); 17.5-19.5" (Lift) Captain: 20.5-22.5" (Flex); 21.5-23.25" (Lift)	SSP: 16.5-18.25" (Flex); 17.5-19.5" (Lift) Captain: 20.5-22.5" (Flex); 21.5-23.25" (Lift)	SSP: 16.5-18.5" (Static TB3) Captain: 20.5-22.5" (Static TB3)	N/A
OPTIONAL POWER ELEVATION	10" Safe Seat Elevation at 3.2 mph	12" iLevel® at 3.5 mph	12" iLevel® at 4.5 mph	12" iLevel® at 3.5 mph	12" iLevel® at 3.5 mph	N/A

- 1) Standard

2) Optional

3) Requires HD Suspension

4) Speed varies with user weight, terrain type
battery charge, battery condition and tire
pressure.

5) Without front riggings

6) Battery capacity may change when adding
power positioning systems.

7) Weight capacity is dictated by the selection
of seating.

8) All measurements +/- 1/8" with 18" x 18" seat,
fixed towers, 300 lbs. weight capacity and solid
tires (add 3/8" with pneumatic tires). Depths
may increase seat-to-floor height. Contact
Customer Service for details. Listed
measurements are in the lowest position;
does not include elevating system unless
specified.

9) 250 lbs. weight capacity with TRU-Balance®

10) A minimum seating size of 20" x 20" is required
with iLevel® on the Q6 Edge® HD; maximum
seating size is 24" x 24"

11) Without foot riggings

12) All Quantum units can be ordered with
specialty controls

13) 6.25 mph speed available only with Q-Logic 3.

14) Battery weight may vary +/- 3 lbs. based on
manufacturer.

15) Reduced to 400 lbs. with iLevel® power
adjustable seat height

16) Seat-to-floor heights are measured with the
unit unloaded. Due to the 4Front's suspension
system, seat-to-floor heights are approximately
0.25" lower when the unit is loaded.

17) 7 mph can only be achieved with optional
4-pole motors and Q-Logic 3. Maximum
speed without 4-pole motors is 6.25 mph
with Q-Logic 3.

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement. The information contained herein is correct at the time of publication. We reserve the right to alter specifications without prior notice.

Quantum® Power Chairs are FDA Class II Medical Devices designed to aid individuals with mobility impairments.



	J4E	J4 HD	J6	4FRONT® 2	4FRONT® 2 HD	Q1450	R-TRAK
HCPCS CODES For more information on HCPCS code coverage, please consult our Product Planning & Reimbursement Center.	Group 2: Single Power (K0835)	Group 2: Single Power HD (K0837)	Group 2: Standard (K0822/K0823), Single Power (K0835)	Group 3: Standard (K0848/K0849), Single Power (K0856/K0857), Multiple Power (K0861)	Group 3: Heavy Duty (K0850), Single Power HD (K0858), Multiple Power HD (K0862)	Group 3: Very HD (K0852/K0853), Single Power Very HD (K0860)	Group 3: Standard (K0848/K0849), Single Power (K0856/K0857), Multiple Power (K0861)
AVAILABLE ON VA CONTRACT	Coming Soon	FSS Contract V797D-70001	FSS Contract V797D-70001	FSS Contract V797D-70001	FSS Contract V797D-70001	FSS Contract V797D-70001	Coming Soon to FSS Contract
WEIGHT CAPACITY ⁷	300 lbs.	400 lbs.	300 lbs.	300 lbs.	450 lbs.	600 lbs.	300 lbs.
MAXIMUM SPEED ⁴	Up to 5 mph	Up to 5 mph	Up to 4 mph	Up to 6.25 mph (3.5 mph when elevated)	Up to 5 mph (2.7 mph when elevated)	Up to 4.5 mph	Up to 7 mph (3.5 mph when elevated)
GROUND CLEARANCE	3" (battery tray)	2.9" (battery tray)	2" (caster arms)	3.14" (battery tray)	3.14" (battery tray)	1.5" (anti-tip wheels)	3.97" (battery tray)
TURNING RADIUS ⁵	20.5"	21"	22"	25.9"	25.9"	26.5"	25.7"
LENGTH ¹¹	35.5"	35.5"	34.74"	39.9"	39.9"	38.25"	34.8"
WIDTH	24.4"	24.5"	23.19"	24"	24"	29.25"	24.2"
DRIVE WHEELS	14" Solid	14" Flat Free	10" Solid Optional: 10" Pneumatic	14" Solid Optional: 14" Pneumatic; 14" Knobby	14" Solid Optional: 14" Pneumatic; 14" Knobby	14" Solid	14" Solid Optional: 14" Pneumatic
CASTER WHEELS	Front: 6" Solid Rear: 6" Solid	Front: 6" Solid Rear: 6" Solid	Front: 5" Solid Rear: 6" Solid	Rear: 9" Flat Free	Rear: 9" Flat Free	Rear: 9" Solid	Front: 9" Solid
ANTI-TIP WHEELS	N/A	N/A	N/A	Front: 4" Flat Free	Front: 4" Flat Free	Front: 4" Solid	Rear: 4" Flat Free
BASE WEIGHT	138.6 lbs.	131.6 lbs.	100.2 lbs.	130 lbs.	130 lbs.	195 lbs.	128 lbs.
BATTERY WEIGHT ¹⁴	38 lbs. each	40 lbs. each	23.4 lbs. each	Group 24: 51.5 lbs. each Group 34: 44.8 lbs. each	Group 24: 51.5 lbs. each	52 lbs. each	51.5 lbs. each
DRIVE WHEEL CONFIGURATION	Mid-wheel Drive	Mid-wheel Drive	Mid-wheel Drive	Front-wheel Drive	Front-wheel Drive	Front-wheel Drive	Rear-wheel Drive
MOTOR PACKAGES	2-pole	4-pole	2-pole	4-pole	4-pole	4-pole	4-pole
BATTERIES ⁶	NF-22 (2)	NF-22 (2)	U-1 (2)	Group 24 (2)	Group 24 (2)	Group 24 (2)	Group 24 (2)
BATTERY CHARGER	8A Off-board	8A Off-board	5A Off-board	8A Off-board	8A Off-board	8A Off-board	8A Off-board
AVAILABLE ELECTRONICS ¹²	120A Q-Logic 3E 75A Q-Logic NE+ (coming soon)	75A Q-Logic NE+	75A Q-Logic NE 75A Q-Logic NE+	120A Q-Logic 3	120A Q-Logic 3	120A Q-Logic 3 90A VR2 6-key	120A Q-Logic 3 75A Q-Logic NE+
SEATING SIZES							
TRU-BALANCE® ⁴	N/A	N/A	N/A	W: 16-22" D: 14-22"	N/A	N/A	N/A
SYNERGY®	N/A	N/A	W: 10-22" D: 10-22"	N/A	N/A	W: 20-32" D: 20-28"	N/A
TRU-BALANCE® ³ STATIC	N/A	N/A	W: 14-20" D: 14-20"	W: 12-22" D: 12-22"	N/A	N/A	W: 12-22" D: 12-22"
TRU-BALANCE® ³ POWER TILT	N/A	N/A	W: 14-20" D: 14-20"	W: 12-22" D: 12-22"	N/A	N/A	W: 12-22" D: 12-22"
TRU-BALANCE® ³ LIFT & TILT	N/A	N/A	N/A	W: 12-22" D: 12-22"	N/A	N/A	W: 12-22" D: 12-22"
TRU-BALANCE® ³ HD STATIC	N/A	N/A	N/A	N/A	W: 20-24" D: 20-24"	N/A	N/A
TRU-BALANCE® ³ HD POWER TILT	N/A	N/A	N/A	N/A	W: 20-24" D: 20-24"	N/A	N/A
TRU-BALANCE® ³ HD LIFT & TILT	N/A	N/A	N/A	N/A	W: 20-24" D: 20-24"	N/A	N/A
TRU-BALANCE® ² HD STATIC	N/A	N/A	N/A	N/A	N/A	W: 20-32" D: 20-28"	N/A
TRU-BALANCE® ² HD POWER TILT	N/A	N/A	N/A	N/A	N/A	W: 20-32" D: 20-28"	N/A
TRU-BALANCE® FLEX SEAT	N/A	N/A	N/A	W: 15-22" D: 15-20"	N/A	N/A	W: 15-22" D: 15-20"
Q4 POWER TILT	W: 16-20" D: 16-22"	W: 20-24" D: 20-24"	N/A	N/A	N/A	N/A	N/A
SEAT-TO-FLOOR HEIGHTS ^{13,16}							
TRU-BALANCE® ⁴ ⁸	N/A	N/A	N/A	N/A	N/A	N/A	N/A
SYNERGY® / STATIC ⁸	N/A	N/A	Synergy®: 15.875 - 17.875" Static: 16.75"	17.75 - 18.75" (Synergy® N/A)	17.75 - 18.75" (Synergy® N/A)	16.5-18.5"	16-18" (Static only)
POWER TILT ⁸	N/A	N/A	16.75"	17.75-19.75"	17.75-19.75"	19.5 - 20.5"	18.5-19.5"
LIFT & TILT ⁸	N/A	N/A	N/A	17.75-19.75"	17.75-19.75"	N/A	17.25-18.25"
TRU-BALANCE® FLEX SEAT ⁸	N/A	N/A	N/A	16.5-18.5"	16.5-18.5"	N/A	16.5-17.5"
Q4 POWER TILT	17.5-19.5"	17.25-19.25"	N/A	N/A	N/A	N/A	N/A
Q4 LIFT AND TILT	17.5-19.5"	N/A	N/A	N/A	N/A	N/A	N/A
Q CAPTAIN SEAT	N/A	SSP: 16.5-18.5" (Flex) Captain: 20.5-22.5" (Flex)	N/A	SSP: 17.5-22.5" (Flex); 17.75-19.75" (Lift) Captain: 21.75-22.75" (Flex); 21.75-23.75" (Lift)	N/A	SSP: 18.5-20.5" (Static TB3) Captain: 22.5-24.5" (Static TB3)	SSP: 15.5-17.5" (Flex); 17.25-18.25" (Lift) Captain: 19.5-21.5" (Flex); 21.25-22.25" (Lift)
OPTIONAL POWER ELEVATION	N/A	N/A	N/A	10" Safe Seat Elevation at 3.5 mph	10" Safe Seat Elevation at 2.5 mph	N/A	10" Safe Seat Elevation at 3.5 mph

quantumrehab.com •     
The information contained herein is correct at the time of publication.
We reserve the right to alter specifications without prior notice.

life.quantumrehab.com
(US) 866-800-2002 • (Canada) 888-570-1113



1 F O R R E H A B P O W E R

Quantum® Power Chairs are FDA Class II Medical Devices
designed to aid individuals with mobility impairments.